

KATALOG 22

PRODUCTS FOR DENTAL PRACTICE
AND LABORATORY



Innovation aus Tradition



since 1949



Innovation from tradition

In turbulent times, we learn again to value virtues that had already seemed to be a bit outdated: reliability, solidarity, perseverance. These are the values that distinguish us as a family-owned company.

Not long after its foundation as a small dental laboratory, Dreve has become an internationally recognized and innovative manufacturer of high-grade products for dental clinics and laboratories. We connect years and years of practical experience with state of the art technical know-how and most recent scientific developments. The result: More than 80 patents, business relations in over 100 countries, and products focusing on you, our customer and most important indicator.

Despite general social distancing, we have managed to move even closer together, and together we work for the successful continuation of the company. One of the many results of our permanent efforts, containing exciting new and innovative products from all areas of the dental world, is in your hands right now.

The fact that we are still able to cope with the tense global situation, including the crisis, keep our optimism and see our chances, is above all owed to you, our customers. At this point we would like to thank you for your loyalty and your trust on behalf of all employees and the Dreve family.



Sincerely yours

Dr. med. dent. Volker Dreve

Guarantee

When manufacturing our products, our staff uses selected materials fulfilling the highest technical requirements. Most careful processing guarantees, in addition to our years of experience, an excellent quality level. Our product systems – which means units, materials and appliances – are perfectly matched, as they are tested and developed in our laboratories every day. Due to the years of experience we are able to guarantee a minimum durability of 2 years from date of production on most of our materials. Of course, we are at your disposal after the acquisition of a product.

Service

A full range of additional services belongs to our performance. Professionally educated dental technicians are at your disposal on our hotline to answer your questions about our products. High-ranking experts from Germany and abroad introduce you to our product systems by means of training courses and seminars. Practical training in our education rooms assumes an important place in the further education programme of Dreve Dentamid GmbH. Information about current presentations can be required day and night on our website www.dreve.com.

Quality






























The quality standard of Dreve products, of course, follows the current level of science. Dreve Dentamid GmbH is certified according to DIN EN ISO 9001, DIN EN ISO 13485 and Directive 93/42/EWG, Annex II. Security for patients, dental technicians and dentists by use of biocompatible products is Dreve’s main focus.

Dreve Dentamid GmbH
Max-Planck-Straße 31
59423 Unna/Germany

Phone +49 2303 8807-37
Fax +49 2303 8807-55
Email dentamid@dreve.de
Internet www.dentamid.dreve.de

Legend

You may already have noticed yourself: We recognize contents of pictures a lot faster than of texts. 60 times faster, to be precise. You will find many new icons in this catalogue, so you can keep facts and numbers in mind easily.

 Hardness, compressive strength, compressive hardness	 Consistency	 Color
 Viscosity	 Film thickness	 Mixing time
 Processing time	 Period inside mouth	 Deflasking after
 Final curing	 Recovery from deformation	 Deformation under pressure
 Shrinkage	 Tensile strength	 Tear resistance
 Elastic modulus	 Flexural strength	 Tensile-, Break-, Bending strength
 Voltage, Operating voltage charger	 Power consumption	 Polymerisation peak
 Pressure, Vacuum	 Temperature range	 Filling capacity
 Dimensions (HxWxD mm)	 Weight	
 Usable blank	 Inner dimensions pressure pot, -chamber (HxWxD mm)	 Time programs, time cycle

NEW: Application tips from our experts




Innovative products from tradition

In perfect shape for you

Impression silicones

Silicone products are one of our core competences. Innovative and economical at the same time: We make this claim to ourselves and our products alike. This is made possible with high-quality raw materials, qualified employees and certified manufacturing processes.

From the transparent impression silicone to the scannable bite registration – our strict quality standards guarantee the consistently high quality that is essential for your daily work.

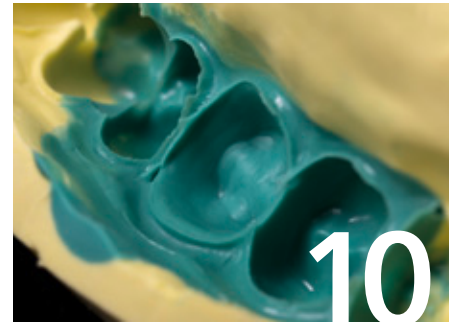

Almost too good to serve as gap filler only

Temporary C&B materials

Fixtemp C&B materials offer outstanding technical values and convince by a setting temperature which is not damaging the dental pulp. They are easy to apply, very precise and reach the soft-elastic phase after just 90 seconds.

In addition to an extension of the color scale, Fixtemp 10:1 now offers a composite that, thanks to improved physical properties, simplifies the fabrication of temporary crowns and bridges as well as inlays and onlays while ensuring superior fit.



Dreve sets standards

3D printing

Leader in detail precision in DLP resin production: With the FotoDent® product world, you benefit from our longstanding expertise with many innovations in medical 3D printing. Perfectly matched systems meet your requirements in CAD/CAM application, enabling you to obtain better results, more precise and longer lasting.

And if you prefer to have it printed: Print@Dreve is your convenient order portal. The dental models can be delivered straight to the laboratory in the highest precision of detail at the push of a button.



Pioneering for decades

Thermoforming technique

We have certainly revolutionised dental technology with the introduction of thermoforming. Still convinced of this process, we continue to work on improving thermoforming technique with high-end equipment and specialised blanks.

Convince yourself of the many possible applications – from individual sports mouthguards to the aligner splint.



Strong alone, stronger together

The resin pouring technique

Functional, aesthetic, custom-fit: With Castdon we offer an unbeatable, complete cold-curing system for the resin pouring technique. In contrast to the hot polymer, dentures made of Castdon resin have a significantly lower shrinkage. In addition, the processing is easier, faster and cleaner.

Are you still unsure in the application? We regularly offer workshops on the topic, held by experienced experts.



Let there be light

Light polymerisation

Light-curing materials are among the most common laboratory helpers in the dental daily laboratory use. Without these and the appropriate equipment many elementary work steps cannot be adequately performed. This is why we constantly research innovations in the field of light curing. The replacement of light tubes by LEDs in our units is a shining example for simplified, more economical work processes.



Runs like a charm

Laboratory silicones

Laboratory silicones have been a fixed component of our product portfolio since the company was founded. Over the years, we have become specialists in dental silicones. We are proud of our product diversification with a wide range of services for everyday work in the dental laboratory – from duplicating to positioner silicone.



The finishing touch

Products and accessories

Small but powerful: Often it is tiny details that distinguish the quality of a piece of work. That's why we dedicate ourselves to the little things with the utmost care and commitment – from the model stump material to the mixing canula. We also set the highest quality standards for our accessories with the guiding principle "Everything from one source". Always with the aim to provide a quality product available for the user.



Dynax® silicones

A whole system for impression taking.

Our Dynax® silicones guarantee both: precision and ease of use. High hydrophilic properties provide excellent impression precision. Our innovative snap-effect enables long processing time but short time in the patient's mouth. The whole impression process takes less than 3 minutes.

- + Snap-effect
- + Optimised hydrophilic properties
- + Scannable surface
- + Improved mechanics



Visual control of the implant impression with Dynax® clear

Dynax® putty

Addition-curing putty impression silicone. Non-sticky, easily kneadable consistency. Suitable for sandwich-, correction double mixing technique and valve border impressions. Combinable with Dynax® light.

2 x 450 ml tubs with dosing spoons	D6101
380 ml cartridge, 5 mixing tips, 1 bayonet-ring	D6150



Dynax® heavy body

Stable, addition-curing "tray-type" silicone. Pressure developing consistency with thixotropic characteristics. Especially suitable for sandwich, correction and double mixing technique. Can be combined with Dynax® light.

2 x 50 ml double cartridges with 6 mixing tips	D6112
8 x 50 ml double cartridges	D6113
380 ml cartridge, 5 mixing tips, 1 bayonet-ring	D6140



Dynax® light

Light bodied, addition-curing correction silicone. Thixotropic, thin flowing correction material (light body) with high detail precision. Suitable for the correction technique can be combined with Dynax® clear, putty, heavy body and mono.

Tip: Dynax® materials are perfectly matched with each other. For best detailed impression we recommend to use Dynax® light only with hydrophilic base material Dynax® putty and Dynax® heavy body.

2 x 50 ml double cartridges with 12 mixing tips	D6103
8 x 50 ml double cartridges	D6104



Impression taking Impression silicones



Impression taking

Dynax® mono

Addition-curing monophase impression silicone. Thixotropic and stable. Suitable for the monophase method, can be combined with Dynax® light.

2 x 50 ml double cartridges with 6 mixing tips	D6106
8 x 50 ml double cartridges	D6107
380 ml cartridge, 5 mixing tips, 1 bayonet-ring	D6120



Dynax® clear

Clear, transparent, addition-curing, monophase impression silicone ensures a perfect control already during the impression-taking. Excelent flowing-characteristics. Especially suitable for the implant impression-taking. Can be combined with opaque correction silicone Dynax® light.

2 x 50 ml double cartridges with 6 mixing tips	D6121
8 x 50 ml double cartridges	D6122



Dynax® putty	1:30 min.	1:30 min.	60 Shore A	kneadable – type 0	sunny yellow
Dynax® heavy body	1:30 min.	1:30 min.	70 Shore A	high – type 1	sky blue
Dynax® light	1:30 min.	1:30 min.	46 Shore A	low – type 3	emerald green
Dynax® mono	1:30 min.	1:30 min.	58 Shore A	medium – type 2	lavender violet
Dynax® clear	1:30 min.	2:00 min.	60 Shore A	medium – type 2	transparent

Impression taking Impression silicones

Zerosil® soft with color indicator

THE alternative to alginates. Now with a chemical color indicator that visually reflects the setting process using the color gradient 1:1. The final color is only achieved when the curing phase has really ended. In addition, processing time and time in the mouth have been reduced by 25 %.

380 ml cartridge	D5118
------------------	-------

Zerosil® soft

Addition-curing impression silicone. High-quality, low-cost alternative to alginate impressions. Impressions can be filled several times. Disinfectable and compatible with all plaster materials. Unlimited storage of impressions. Especially suitable for the situation impression.

2 x 500 ml tubes	D48411
2 x 875 ml tubes	D4841



Your silicone for any application

	Dynax® clear	Dynax® heavy body	Dynax® light	Dynax® mono	Dynax® putty	Zerosil® soft
Indication						
Crown and bridge impressions						
Double mixing technique	•	•	•	•	•	
Correction	•	•	•	•	•	
Sandwich	•	•	•	•	•	
Monophase	•	•	•	•	•	
Inlay and onlay impressions						
Double mixing technique	•	•	•	•	•	
Correction	•	•	•	•	•	
Sandwich	•	•	•	•	•	
Monophase	•	•	•	•	•	
Pick-up impressions e. g. for implants						
Correction	•	•		•		
Monophase	•	•		•		
Functional impressions						
Function				•		
Fixation impressions						
Monophase		•		•		
Relining impressions						
Monophase			•	•		
Pressure spot indicator						
Monophase			•	•		
Impressions for temporaries						
Monophase				•		•
Situation				•		•
Situation impressions						
Monophase						•
Situation						•
Impressions of the opposite jaw						
Monophase						•



Situation impression with Zerosil® soft with color indicator



Silicones need more than 30 minutes to reach their perfect characteristics. Please make sure that you do not fill the silicone impression immediately after taking it out of the patient's mouth.

Impression taking Bite registration

StoneBite®

Addition-curing silicone with an extreme final hardness of 48 Shore D. Absolutely precise and reliable. Especially suitable for bite registrations and registration keys for support pins.

2 x 50 ml double cartridges, 12 mixing tips, 12 application tips	D504872
4 x 50 ml double cartridges	D504873

StoneBite® scan

Scannable bite registration material based on A-silicone. High final hardness, thixotropic and stable. Very easy to cut and mill. Specially developed for use in the CAD/CAM area.

2 x 50 ml double cartridges, 12 mixing tips, 12 application tips	D504870
------------------------------------------------------------------	---------

Regofix transparent


Transparent bite registration material, free from filling materials, not opaque to x-rays. Easy to cut with a high final hardness and spontaneous recovery from deformation. Suitable as a matrix material for light-curing anterior and posterior tooth composites and for guided bite registration.

2 x 50 ml double cartridges, 12 mixing tips	D504876
8 x 50 ml double cartridges	D504877

Regofix opak

Bite registration material with very good cutting and milling characteristics. Fast, hard, precise. Suitable for all kinds of bite registrations and registration keys for support pins.

2 x 50 ml double cartridges, 12 mixing tips	D504874
8 x 50 ml double cartridges	D504875

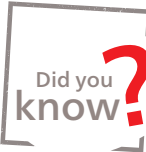

Shore D is the measuring unit for the determination of the hardness of non-flexible materials, if beyond 90 Shore A.



StoneBite®: Exact delineation, precise processing, long-term stability



Zerosil® soft with color indicator	1:00 min.	2:00 min.	40 Shore A	medium – type 2	plum to orange
Zerosil® soft	1:30 min.	1:30 min.	40 Shore A	medium – type 2	orange
StoneBite®	0:30 min.	1:20 min.	48 Shore D	medium – type 2	orange
StoneBite® scan	0:30 min.	1:20 min.	32 Shore D	medium – type 2	grey
Regofix transparent	0:30 min.	1:30 min.	70 Shore A	medium – type 2	transparent
Regofix opak	0:30 min.	1:30 min.	90 Shore A	high – type 1	blue

Impression taking Accessories

Mixing tips

Mixing ratio 1:1
Inner-Ø 4.2 mm
Color mint green



Dynax® light

40 pieces	D4970X7
100 pieces	D49701X7

Mixing ratio 1:1
Inner-Ø 5.4 mm
Color pink



StoneBite®, StoneBite® scan

40 pieces	D4970X1
100 pieces	D49701X1

Mixing ratio 1:1
Inner-Ø 4.2 mm
Color yellow



Dynax® light

40 pieces	D4970X5
100 pieces	D49701X5

Tips

For mixing tips for intra-oral application.

transparent (e. g. for Dynax® clear)
yellow (e. g. for Dynax® light)



40 pieces, transparent	D49602
40 pieces, yellow	D49604

Mixing ratio 1:1
Inner-Ø 6.5 mm
Color green



Dynax® clear, Dynax® heavy body
Dynax® mono, Regofix

40 pieces	D4970X2
100 pieces	D49701X2

Mixing ratio 5:1
Inner-Ø 6.0 mm
Color yellow



Zerosil® soft, Dynax® putty
Dynax® heavy body, Dynax® mono

15 pieces	D4982
50 pieces	D4981

Application tips

For mixing tips for intra-oral application.

transparent (e. g. for StoneBite®)



40 pieces, transparent	D49603
------------------------	--------

Bayonet-ring

Attachment for the use of 380 ml cartridges.



Bayonet-ring, yellow	D4980
----------------------	-------

DS 50 Injector 1:1

For applications of dental laboratory silicones for S 50 double cartridges.

Mixing Ratio 1:1

DS 50 Injector 1:1	D5007
--------------------	-------



Temporary C&B materials

Fixtemp C&B 4:1

Automatically mixable composite in a mixing ratio of 4:1 for the manufacturing of temporary crowns, bridges, inlays and onlays. High breaking strength, no temperature development which damages the pulp at intraoral setting. Fast and secure application through safety-cartridge.

Fixtemp C&B 10:1

Automatically mixable composite in a mixing ratio 10:1 for the preparation of temporary crowns and bridges as well as inlays and onlays. Highest fitting accuracy with the least abrasion. Improved physical properties for easier grinding and polishing.

Color	Fixtemp C&B 4:1		Fixtemp C&B 10:1	
	50 ml dc + 10 mt	5 x 50 ml dc + 10 mt	50 ml dc + 10 mt	5 x 50 ml dc + 10 mt
A1	D50391		D60391	
A2	D50392	D503925	D60392	D603925
A3	D50393	D503935	D60393	D603935
A3,5	D50394		D60394	
B1	D50395		D60395	
Bleach X	D50396		D60396	
D2			D60399	



Do not miss the point of time to take Fixtemp C&B out of the patient's mouth when still in the flexible phase. Every resin has a shrinkage!



Convincing results with Fixtemp C&B

Fixtemp C&B 4:1	0:45 min.	~ 2:00–3:00 min.	~ 6:00 min.	≥ 200 MPa	≥ 60 MPa	–
Fixtemp C&B 10:1	0:45 min.	~ 2:00–3:00 min.	~ 6:00 min.	≥ 300 MPa	≥ 80 MPa	–

Mixing tips

Mixing ratio 4:1/10:1
Inner-Ø 3.2 mm
Color light blue



Fixtemp C&B 4:1
Fixtemp C&B 10:1

40 pieces	D4970X4
100 pieces	D49701X4

DS 50 Injector 4:1/10:1

For applications of dental provisional C&B materials from S 50 double cartridges.

Mixing ratio 4:1/10:1

DS 50 Injector 4:1/10:1	D5008
-------------------------	-------



FotoDent® biobased model

Our first green 3D printing resin

Sustainability, climate, environment are as topical issues as ever and will remain so. We take this development into account and can now present the first impressive result to the dental world with the newest member of the FotoDent® family:

FotoDent® biobased model is a light-curing resin that consists of 50 percent renewable raw materials. The ecological advantage of bio-based raw materials over petroleum-based raw materials is that they have less impact on the climate during the production of the resin: If the carbon is bound from biomass, it does not turn into carbon dioxide. This is where plants convert carbon dioxide into a valuable raw material. This reverses the trend that humanity is increasingly emitting more carbon than carbon dioxide into the atmosphere.

Did you know?

The renewable raw materials are parts of wheat and maize that are not used for food production. The production of these herbal ingredients does not compete with food production. They are also genetically unmodified.

- ⊕ Reduced impact on the atmosphere
- ⊕ Free from carcinogenic substances
- ⊕ Shorter cleaning times



FotoDent® product world

Generative production from one source.

Already for decades Dreve Dentamid GmbH with their digital dental products stand for competence and innovation in the field of medical 3D printing. Just like it is with the classical portfolio, the system concept plays a decisive part: resins and post-curing units, all from one source, are perfectly matched. Thus we provide for specific processes in the dental laboratory.

No matter whether it is lab, prosthetics or orthodontics: only products complying with the high quality stipulations of dental applications and medical devices during the entire production chain, will leave our factory. Results from the FotoDent® product world speak for themselves and add to the success of every future-oriented dental business.



Biobased dental models

3D printing



Dental drill guides




Dental gingiva masks



Castable moldings



Indirect bonding trays

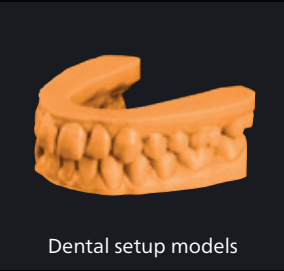

Dental setup models



Denture bases



Dental models



Dental individual trays



Don't forget!
Shake the bottle well before you pour the resin into the vat.

3D printing FotoDent® resins


FotoDent® denture 385 nm
Light-curing resin for the production of denture bases. The material provides for accurate results with an outstanding dimensional stability. FotoDent® denture is MMA-free and stable in color.

pink transparent, 1.0 kg bottle	D35500
pink opaque, 1.0 kg bottle	D35501



FotoDent® IBT 385 nm
Light-curing resin for the manufacturing of transparent indirect bonding trays. Monitoring of correct bracket positions is possible at any time – before and after fixing. Simple removing due to optimal flexibility. Usable with all common bracket systems.

1.0 kg bottle	D35120
---------------	--------

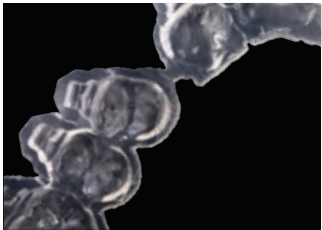

FotoDent® guide 385 / 405 nm
Light-curing resin for the manufacturing of dental drill guides. Drill guides made from FotoDent® guide are dimensionally stable and biocompatible.

385 nm, clear transparent, 1.0 kg bottle	D35650
405 nm, blue transparent, 1.0 kg bottle	D35600



FotoDent® tray 385 / 405 nm
Biocompatible light-curing resin for the manufacturing of dental impression trays. High reactivity at low viscosity. The finished impression trays have a homogeneous surface, which significantly reduces manual reworking.

1.0 kg bottle	D35300
---------------	--------



FotoDent® cast 385 / 405 nm
Light-curing resin for the production of burn-out molded parts for the dental casting technique. The material is highly compatible: it can be used with all commercially available embedding materials.

1.0 kg bottle	D35100
---------------	--------



FotoDent® setup 385 / 405 nm
Light-curing resin for the manufacturing of orthodontic set-up models. The material is above all characterised by its increased building speed. The set-up models are very well suitable for the thermoforming technique.

1.0 kg bottle	D35900
---------------	--------



3D printing FotoDent® resins


FotoDent® biobased model 385 nm
Light-curing resin that consists of 50 percent renewable raw materials. The material properties are of course ideal for model production (also in Aligner workflow). FotoDent® biobased model can be cleaned with FotoClean.

1.0 kg bottle	D35450
---------------	--------



FotoDent® model2* 385 nm
Especially developed for new generations of printers with 385 nm: Light-curing resin for the manufacturing of dental models for orthodontics and prosthetics. Faster printing process due to higher output, still with the same precision.

beige opaque, 1.0 kg bottle	D354002
black, 1.0 kg bottle	D354003
white, 1.0 kg bottle	D354004

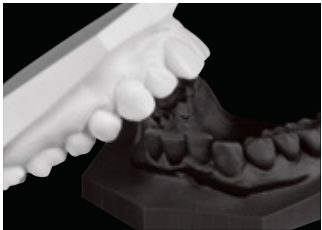

FotoDent® gingiva 385 nm
Light-curing resin for the manufacturing of dental gingiva masks. Permanently soft and flexible. Perfectly combinable with working models made of FotoDent® biobased model and model2 resins.

1.0 kg bottle	D35850
---------------	--------



FotoDent® denture	> 2.000 MPa	> 80 MPa	–	pink transparent, pink opaque	0.4–0.6 Pa s
FotoDent® IBT	–	–	≥ 45 %	clear transparent	0.7 ± 0.2 Pa s
FotoDent® guide 385 nm	≥ 1,900 MPa	≥ 90 MPa	≥ 9 %	clear transparent	0.7 ± 0.2 Pa s
FotoDent® guide 405 nm	≥ 1,700 MPa	≥ 75 MPa	≥ 10–15 %	blue transparent	0.65–1.05 Pa s
FotoDent® tray	≥ 2,000 MPa	≥ 75 MPa	> 5 %	green	~ 1 Pa s
FotoDent® cast	≥ 2,000 MPa	≥ 100 MPa	7.5–11 %	red transparent	< 0.3 Pa s
FotoDent® setup	≥ 2,000 MPa	≥ 95 MPa	≥ 5 %	maize yellow	1.0–1.5 Pa s
FotoDent® biobased model	≥ 2,200 MPa	≥ 95 MPa	≥ 5%	beige opaque	~ 0.6 Pa s
FotoDent® model2	≥ 1,900 MPa	≥ 85 MPa	≥ 8%	beige opaque	1.0 ± 0.2 Pa s
FotoDent® model2	≥ 2,000 MPa	≥ 100 MPa	≥ 10%	black	0.75 Pa s
FotoDent® model2	≥ 1,950 MPa	≥ 85 MPa	≥ 8.5%	white	0.7 ± 0.2 Pa s
FotoDent® gingiva	–	–	> 40 %	pink	~ 2 Pa s

*Please note: Depending on the selected color, the printing times of the models may differ.

3D printing Cleaning

Bandelin RK 100 H

High quality, robust ultrasonic cleaning unit with stainless steel tub. The compact and easy-to-use device has a timer from 1 to 15 minutes or runs continuously. Heating function for the temperature range between 30 and 80°C.



Inner dimensions:
H 100 mm
W 240 mm
D 140 mm

Bandelin RK 100 H	82100
-------------------	-------


- ⊕ Ultrasonic power: 320 W
- ⊕ Infinitely adjustable temperature
- ⊕ Recommended for use with FotoClean



Insert basket made of stainless steel

Dimensions: H 52 x W 175 x D 70 mm.

1 piece	82103
---------	-------

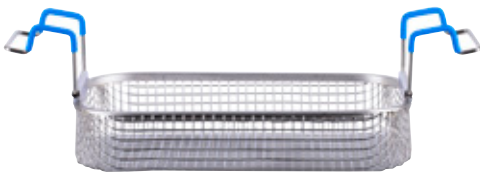

Stainless steel perforated lid

Suitable for 2 beakers.

1 piece	82101
---------	-------

Beaker

600 ml, Ø 85 mm.

1 piece	82102
---------	-------



Bandelin RK 100 H	230 V / 50 Hz
	2.0 l

3D printing Cleaning / Curing

FotoClean

FotoClean is a ready-to-use, water-based cleaning mixture for use in the heated ultrasonic cleaning device Bandelin for cleaning additively manufactured components.

2 x 1 l bottle	50000
5 l canister	50001



Inner dimensions
polymerisation chamber:
H 65 mm
W 150 mm
D 150 mm



PCU LED N₂

Laboratory unit on LED-basis for the curing of 3D printed parts. Ensuring of mechanics and biocompatibility through effective deep curing. Two different curing environments selectable depending on the application area. The additionally eligible nitrogen environment of PCU LED N₂ ensures cured components without inhibition layer – for laboratory and medical devices.

PCU LED N ₂	4317
------------------------	------

- ⊕ Open system with 10 programmable memory spaces
- ⊕ Simple operating concept with electronic control system
- ⊕ Logging and monitoring of process parameters
- ⊕ Selectable curing environments

PCU LED N ₂	100–240 V / 50–60 Hz, 0.7 A	9.3 kg	110x389x276 mm
			red

Thermoforming technique by Dreve

Pioneering for decades.

We have certainly revolutionised dental technology with the introduction of thermoforming. Still convinced of this process, we continue to work on improving thermoforming technique with high-end equipment and specialised blanks.

Convince yourself of the many possible applications – from individual sports mouthguards to aligners.



Our first mouthguard with a 3D-printed insert for best protection and fit.
Recommended by world boxing champion Christina Hammer.


Thermoforming blanks

All Dreve thermoforming blanks are round as standard. The article numbers listed in the catalog refer to **round blanks with a diameter of 120 mm**.

Do you need round blanks with a diameter of 125 mm?

Simply add an “S” after the article number.

Example: D420001 (Biolon 120 mm) / D420001S (Biolon 125 mm)

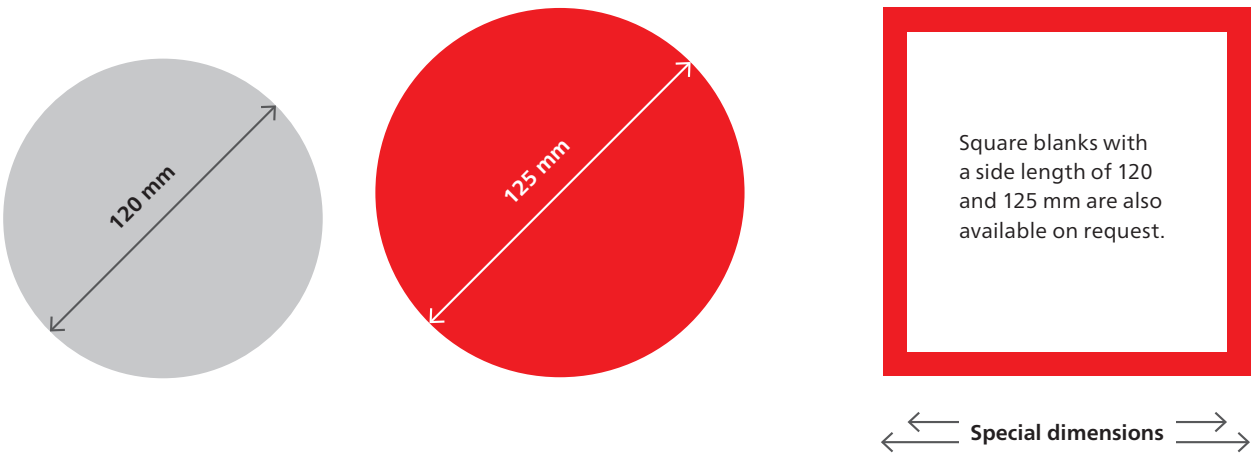

Don't forget!
Pay attention to the recommended foil diameter of your thermoforming unit!

Tip from the pro
Always pre-heat your unit for at least one minute without thermoforming foil.

Application	Product	Strength				
Aligner	Biolon	0.5 mm	0.625 mm	0.75 mm		
Base plates	Drufoplast	2 mm				
	Drufoplast incl. tray handle	3 mm				
Bleaching	BioBleach® hard	0.5 mm				
	BioBleach® soft	1 mm				
	Drufosoft®	1 mm	1.5 mm			
Matrix	Drufolen H	0.6 mm				
	Drufolen W	0.45 mm	0.65 mm	0.8 mm		
Custom-made trays	Drufoplast incl. tray handle	3 mm				
Orthodontic devices	Biolon	2 mm	3 mm			
Mouthguards	Biolon	2 mm				
	Drufosoft® transparent	2 mm	3 mm	4 mm		
	Drufosoft® colour	3 mm				
Space maintainers	Drufolen W	0.65 mm	0.8 mm			
Positioner / Retainer	Biolon	0.75 mm				
	Drufosoft® transparent	1.5 mm	2 mm	3 mm	4 mm	5 mm
	Drufosoft® colour	3 mm	4 mm			
	Drufosoft® pro	3 mm				
Production of temporaries	Biolon	0.5 mm	0.75 mm			
	Drufolen H	0.6 mm	0.8 mm			
X-ray splints	Biolon	1 mm	1.5 mm	2 mm		
	Drufoplast incl. tray handle	3 mm				
Splint therapy	Biolon	0.75 mm	1 mm	1.5 mm	2 mm	3 mm
	Drufosoft® pro	3 mm				
	Kombiplast	2 mm	3 mm	4 mm		

Thermoforming technique Blanks

BioBleach® hard

Thermoforming blank made of polyethylen especially developed for bleaching and medical trays. Slightly flexible, also suitable for situations with slight undercut areas. Delivery including insulating foil for high, clear transparency.

$\frac{3}{8}$ mm	pieces	item No.
0.5	25	D420015

BioBleach® soft

Transparent thermoforming blank made of vinyl acetate copolymer with a hardness of 95 Shore A. For the production of soft, bleaching/medication splints extending the gingival margin.

$\frac{3}{8}$ mm	pieces	item No.
1	20	D4300

Biolon

The thermoforming blank made of PETG (polyethylene terephthalate glycol) is supplied with an insulating film for the purpose of high transparency. Since it is tough-elastic and can be built up with acrylate, it is particularly suitable for aligner, bite and occlusion splints, as well as for orthodontic appliances.

$\frac{3}{8}$ mm	pieces	item No.
0.5	25	D420001
	50	D420002
0.625	20	D420019
	40	D420020
0.75	20	D420003
	40	D420004
1	25	D420005
	50	D420006
1.5	30	D420007
	60	D420008
2	20	D420009
	60	D420010
3	10	D420011
	50	D420012

Bleaching/Medicament splint made of BioBleach® hard. Cenit L used as a spacer



Bleaching/Medicament splint made of BioBleach® soft. Drufolen W used as a spacer



Aligner made of Biolon



Thermoforming technique Blanks



Biolon for Vacfomat V9 2.0

Tough elastic hard plastic out of PETG (Polyethylene terephthalate glycol). Up to 9 aligners can be produced per blank in one thermoforming process. Biolon is equipped with an insulating foil to give the aligner a highly transparent appearance.

$\frac{3}{8}$ mm	Ø mm	color	pieces	item No.
0.5	354	clear transparent	50	D420021
0.625				D420022
0.75				D420023
1				D420024



9 blanks in one thermoforming process: Biolon for Vacfomat V9 2.0

Drufolen H

Thermoforming blank made of polyethylene. Due to its hard-flexible consistency and residue-free combustion it is especially suitable for the production of molding parts, transfer caps and stencils for the production of temporaries.

$\frac{3}{8}$ mm	Ø mm	color	pieces	item No.
0.6	42	natural	150	D4237442
0.6	120	natural	40	D423641
0.8	120	natural	40	D423681

Template for the production of provisoricals made of Drufolen H



Drufolen W

Thermoforming blank made of polypropylene. Due to its soft-flexible consistency and residue-free combustion it is especially suitable for the production of molding parts and matrix with undercut areas and space maintainers used for medical splints.

$\frac{3}{8}$ mm	Ø mm	color	pieces	item No.
0.45	42	natural	150	D42384N42
	120	natural	40	D423441N
0.65	42	natural	100	D42386N42
	120	natural	40	D423461N
	120	blue	40	D423461B
	42	natural	100	D42388N42
0.8	120	natural	40	D423481N
	120	blue	40	D423481B
	120	natural	30	D423401N

Bleaching/Medicament splint with Drufolen W used as a spacer



Sample selection thermoforming blanks

The complete range of the Dreve thermoforming blank program contains all the information you need about materials, layer thicknesses and accessories. Packed in a handy binder.

Binder	D3272
--------	-------



Blank trays

Compact overview of the complete range of blanks including item numbers, colors as well as heating and cooling times.

Sample booklet	D770178
----------------	---------



Thermoforming technique Blanks

Drufosoft®

Thermoforming blank made of ethyl vinyl acetate. Soft elastic with high recovery from deformation with a final hardness of 82 Shore A. Can be laminated and is ideally suited for sports mouthguards, positioners and duplication forms.

Please note:
Drufosoft® transparent and colour are available with special dimensions on request.

Drufosoft® colour



mm	color	pieces	item No.
1	transparent	20	D4245
		50	D42451
1.5		20	D4252
		50	D42531
2		15	D4250
		50	D42511
		10	D4248
		30	D42481
neon red		3	D42483X2
		10	D4248X2
neon green		3	D42483X3
		10	D4248X3
neon yellow		3	D42483X4
		10	D4248X4
blue		3	D42483X5
		10	D4248X5
yellow	3	D42483X6	
	10	D4248X6	
black	3	D42483X7	
	10	D4248X7	
white	3	D42483X8	
	10	D4248X8	
pink	3	D42483X9	
	10	D4248X9	
3	red	3	D42483X17
		10	D4248X17
turquoise	3	D42483X18	
	10	D4248X18	
gold	3	D42483X19	
	10	D4248X19	
silver	3	D42483X20	
	10	D4248X20	
blue opaque	3	D42483X21	
	10	D4248X21	
green	3	D42483X22	
	10	D4248X22	
purple	3	D42483X23	
	10	D4248X23	
light blue	3	D42483X24	
	10	D4248X24	
bordeaux	3	D42483X25	
	10	D4248X25	
orange	3	D42483X26	
	10	D4248X26	
4	transparent	15	D42472
		75	D424901
5		15	D4225

Thermoforming technique Blanks

Drufosoft® pro

Transparent thermoforming blank made of ethyl vinyl acetate. Hard-flexible occlusion splints with high recovery from deformation with a final hardness of 95 Shore A. Can be laminated and is perfectly suitable for occlusion splints, positioners and bruxism splints.

mm	pieces	item No.
3	10	D4230

Occlusal splint made of Drufosoft® pro




Drufoplast

Red-opaque blanks made of polystyrene. Extremely break-resistant. Especially suitable for base plates for bite-takings and set-ups.

mm	pieces	item No.
2	15	D4210Z

Base plate with Drufoplast

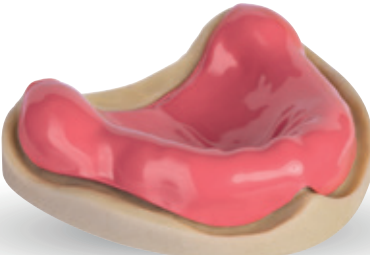

Drufoplast incl. tray handle

Transparent thermoforming blank made of polystyrene extremely high break-resistance. Especially suited for custom-made impression trays and x-ray splints.

mm	pieces	item No.
3	10	D4214
3	100	D4215

Transparent individual tray made of Drufoplast



Kombiplast

Highly transparent two layer thermoforming blank made of PETG and EVA. High wearing comfort of occlusion splints due to hard surface and soft inner side. Bonds with acrylate for the individual adjustment.

mm	pieces	item No.
2 (soft 1 : hard 1)	10	D420016
3 (soft 1.5 : hard 1.5)	10	D420017
4 (soft 1 : hard 3)	10	D420018

Occlusal splint made of Kombiplast



Thermoforming technique Mouthguard

Dreve Mouthguard

A must in every sports bag. The professional mouthguard with optimum fitting and best protection. No impact on talking, breathing and swallowing. Select the appropriate mouthguard according to the sports discipline: By using different materials and material thicknesses, the mouthguard is customised according to the sports activity. X-ray opaque stripes stand for safety. So does the Dreve Mouthguard logo, ensuring that the mouthguard was produced in a certified dental lab. Our sales channel through certified dentists and dental technicians trained by Dreve stands for constantly high quality unrivalled for many years.

- ⊕ Customised fitting by the dentist
- ⊕ Individually manufactured in the dental lab
- ⊕ Optimum shock absorption and firm fitting in the mouth



Dreve Junior

Sports mouthguard for children and adolescents up to the age of 14.

- 1 x Drufosoft® transparent (2 mm)
- 1 x Drufosoft® colour (3 mm)

Dreve Senior

Sports mouthguard for adults, for most hobby sports disciplines like Handball, Basketball, Squash, American Football etc.

- 1 x Drufosoft® transparent (3 mm)
- 1 x Drufosoft® colour (3 mm)



Dreve Elite

Sports mouthguards for martial arts (with soft inlay) e.g. Boxing, Judo, Wrestling, Rugby.

- 2 x Drufosoft® transparent (3 mm)
- 1 x Drufosoft® colour (3 mm)



Dreve Mouthguard Professional 3D

Sports mouthguard with a 3D printed insert.

- 1 x Drufosoft® transparent (3 mm)
- 1 x 3D Protector inlay
- 1 x Drufosoft® colour (3 mm)



Thermoforming technique Mouthguard

Dreve Mouthguard professional 3D

A miracle of bionics: In close exchange with sports scientists and professional athletes, we have succeeded in combining two different methods – traditional thermoforming technology and digital printing technology – into one product that exceeds the previously available mouthguards in terms of protection and comfort.



Dreve Mouthguard professional 3D Starter Set

Complete set for the production of two “next level” mouthguards with a 3D-printed, flexible insert for the best protection and accuracy of fit.

Included:

- 4 x Drufosoft® 3.0 mm
- 4 x 3D Protector inlays (1 x / type)
- 1 x 20 ml Protector Adhesive
- 1 x 20 ml Finishing Liquid
- 1 x Mouthguard Box



All items are also available separately.


Dreve Mouthguard professional 3D Starter Set	D44600
----------------------------------------------	--------



3D Protector inlays

3D-printed inlays (8 pieces) for the production of professional 3D mouthguards. Choose between Basic (0.6 and 1 mm) and Premium (0.9–1.5) depending on the required protection.

mm	Size	Pieces	Artikelnr.
Basic	L	8	D44601
Basic	XL	8	D44602
Premium	L	8	D44603
Premium	XL	8	D44604

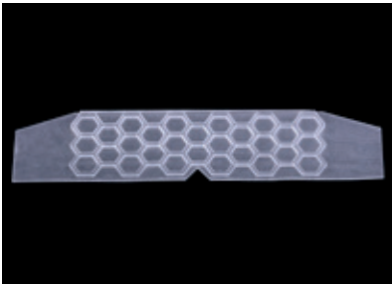

The honeycomb structure of our new 3D Protector inlays is inspired by bionics. Bionics describes the science of solving technical problems based on the example of biological functions.

Protector Adhäsiv

For connecting the first Drufosoft® blank with a 3D Protector insert.

20 ml bottle	D44605
--------------	--------



Thermoforming technique Mouthguard

Dreve Mouthguard Box

For the hygienic and safe storage of mouthguards.

5 pieces, red	D4444
5 pieces, black with imprint	D44606



Would you like to produce your “next level” mouthguard?
Find the detailed, illustrated step-by-step instructions on
the-next-level-mouthguard.com

X-ray opaque stripes

For insertion into Dreve Mouthguards. Allows the detection of the mouthguard on X-ray images, for example in unconscious patients.

15 pieces	D4445
-----------	-------



Dreve Mouthguard Label

Insertion for Dreve Mouthguards. Issued exclusively to certified laboratories. Proof of outstanding quality regarding production and efficiency of the mouthguard.

1 sheet/60 pieces	D44500
-------------------	--------



Thermoforming technique Units

Vacfomat V9 2.0

Vacuum-based thermoforming unit, specially developed for the production of aligner splints. Up to 9 aligners can be thermoformed at the same time in just one operation. Easiest operation as well as constant digital and acoustic control of all work phases ensure problem-free work even with little prior knowledge.

Vacfomat V9 2.0	D3282
-----------------	-------



- ⊕ Vertical molding
- ⊕ Powerful vacuum pump
- ⊕ Up to 108 aligner blanks / hour



Reliably precise results in industrial quality with the Vacfomat V9 2.0








Magnetic pliers V9 2.0

Gripping pliers for safe removal and integration of the model tray in the Vacfomat V9 2.0.



The model tray can be safely removed from the device and reintegrated without direct contact thanks to the magnetic locking of the gripping pliers.

					
Vacfomat V9 2.0	220–240 V / 50 Hz	2,000 W	-0.8 bar	45 kg	415x488x601 mm

Thermoforming technique Units

Druformat scan


Very easy to handle high-end thermoforming unit with a maximum in user-friendliness and consistently precise final results. Built for all applications in dental thermoforming technique. Easy to use keypad, integrated bar code scanner identifies blanks and a working pressure of up to 6 bar. Additional safety factors for professional results. Two-hand operation and LEDs that light up when heating is in operation.

Content:

- Base range of blanks
- Drain pot for granules
- Air connection hose (1 m)

red	D3300
silver	D3300X1
Special color from RAL-Color-Chart	D3300X2


from diameter 120 mm

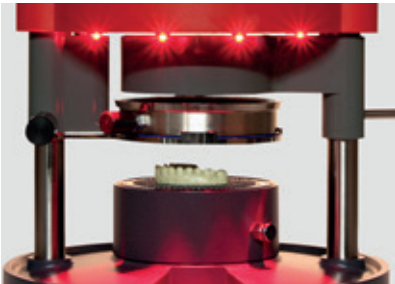

Aligners, mouthguards, bruxism splints: The thermoforming allrounder Druformat scan belongs in every professional lab



- ⊕ Perfect combination of technology, software and design
- ⊕ Very easy operation
- ⊕ High working-pressure for all applications
- ⊕ Multi-functional professional unit



Pressure thermoforming machines work with a forming pressure which is up to 12 times as high as the forming power of vacuum units.

Thermoforming technique Units

Drufosmart scan

Pressure forming unit with integrated bar code scanner to identify blank codes and a working pressure of up to 4 bar. Display with self-explanatory control panel, vertical blank impression for best results. Due to the 2-hand operation and automatic switch-off function of the heating element safe working conditions are ensured. Covers the entire scope of thermoforming technique.

Content:

- Base range of blanks
- Air connection hose (1 m)

red	D3400
-----	-------


from diameter 120 mm



Drufosmart

Basic pressure thermoforming unit with 2.5 bar working pressure. The better alternative in comparison to vacuum units with about three times higher forming pressure. Semi-automatic operation for easy adjustment of the corresponding heating times 2-hand operation and automatic switch off function of the heating element ensure safe working conditions.





Content:

- Base range of blanks
- Air connection hose (1 m)

red	220–240 V / 50 Hz	D3200
red	115 V / 60 Hz	D3200A


from diameter 120 mm



					
Druformat scan	115–240 V / 50–60 Hz	max. 315 W	2–6 bar (19–87 psi)	16 kg	490x290x290 mm
Drufosmart scan	115–240 V / 50–60 Hz	max. 275 W	2–4 bar (19–58 psi)	15 kg	420x320x290 mm
Drufosmart	220–240 V / 50 Hz 115 V / 60 Hz	max. 275 W	2.5 bar (36 psi)	15 kg	420x320x290 mm

Workout set thermoforming technique


Set of grinding tools for the workout and finishing of all thermoformed parts, assorted in a robust and practical table stand, made of stainless steel.

Content:

- 1 x square cutter* (D4374V)
- 1 x cutter* (D4374S)
- 1 x UltraTrimm (medium)* (D43752)
- 1 x UltraTrimm (light)* (D43753)
- 1 x grinding cap holder (size 2)* incl. grinding cap (D38230T2/D38230S2)
- 1 x grinding cap holder (size 1)* incl. grinding cap (D4458/D4461)
- 2 x mandrels incl. support discs
- 1 x hard rubber, flame shaped
- 1 x brush made of goat hair
- 1 x cotton buff

Items marked with * are separately available.

Workout set	D3273
-------------	-------

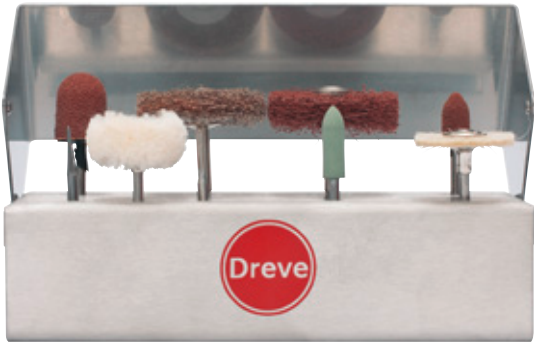

Foil scissors

Stable special scissors for use in the thermoforming technique, with ridged and bent cut.

Foil scissors	D4369
---------------	-------



Cutters

Carbide burs for quick and clean separating and cutting or finishing of all hard thermoforming materials.

Square cutter without toothing	D4374V
Piercing cutter with spiral cut	D4374S



Finishing works with the workout set thermoforming technique



Square cutter Piercing cutter

Finishing Liquid

Polishing fluid for removing roughness from Drufosoft® materials.

250 ml bottle	D4466N
---------------	--------



Fillin

Reusable block-out material on silicone basis, for covering undercut areas on the model.

60 g strip, nature	D42541
--------------------	--------



Accurate work with Fillin

Silicone grease

Special care agent for lubricating rubber seals, for example o-rings on thermoforming units and pressure pots.

35 g tube	D3350
-----------	-------



You live and learn: Dreve Mouthguard workshop

The range of topics of our workshops offered is as diverse as the dental work routine. Thanks to this broad range, dentists, dental technicians, trainees and students benefit from our know-how and numerous practical tips.

With the hands-on courses, you will have the opportunity to learn our range of systems in a well-founded and practice-oriented way. Great emphasis is on “learning by doing”. That is why we attach great importance to practical work in our workshops. That way, you can use what you have learned in the workshop immediately in your daily work and be more efficient in the long run.

Information / Registration:

dentamid.dreve.de/service/workshops



MaxiFresh™

Cleaning tablet for the outer cleaning of dental splints. Used daily it guarantees perfect hygiene and reduction of decoloration. Does not contain any abrasives and thus extends considerably the life span of splints and sports mouthguards.

24 tablets with bowl	D4100
24 tablets refill set	D4101



MaxiFresh™ with bowl

Granules

For embedding unused model parts during the thermoforming process.


Round granules; 1.2 kg bottle	D4282
Steel granules; 1.2 kg bottle	D4283
Refined steel granules; 1.2 kg bottle	D4284



Funky tool

Funky tool is an accessory unit for Drufoformat and Drufosmart scan thermoforming units for an individual color design of Drufosoft® blanks.

Case with different inserts	D3229
-----------------------------	-------

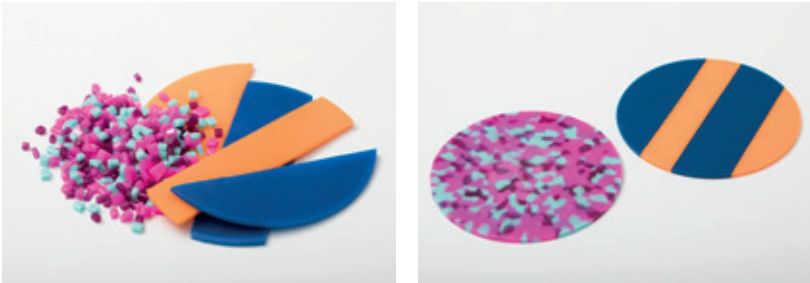

With the appropriate Funky tool insert whole teams can be equipped with mouthguards in club colors

Pressure polymerisation insert

This insert allows an extra use of Drufoformat units as pressure pot for cast flasks. The pressure polymerisation insert is to be placed onto the base plate of the Drufoformat. The desired pressure can be adjusted individually between 2 and 6 bar.

Pressure polymerisation insert	D3227
--------------------------------	-------



Denture production with Castdon resin

Functional, aesthetic, custom-fit.

With Castdon we offer an unbeatable, complete cold-curing system for the resin pouring technique. In contrast to the hot polymer, dentures made of Castdon resin have a significantly lower shrinkage. In addition, the processing is easier, faster and cleaner.



Results will be even more precise if resin is cured under pressure.



Align

- Maximum dorsal model height
- Create wax modeling as usual
- Functional ridges must be free of wax
- Align model in flask



Fix

- Fix model with wax on the base plate
- Lock the flask top with clamps
- Check correct closure



Mix

- If necessary, shake Silifform before use
- Measure 150 ml for each component and mix homogeneously with spatula



Fill

- Fill flask completely from above in a thin stream



Curing

- Ideally, cure the silicone under the same pressure as during the polymerisation of the casting resin (e. g. Polymax)
- Curing about 30 minutes



Demoulding

- Loosen the bottom plate and demould the model



Cut out

- Two sprues are always stung in the dorsal area of the alveolar ridge
- Create a channel in the upper jaw centrally to the base



Repositioning

- Remove wax from teeth with stainless steel bottom plate
- Roughen teeth and attach retention
- Reposition teeth in the duplicating form
- The moist surface of Silifform does not require bonding of the teeth



Bonding

- Wetting with Castdon Bonding strengthens the bond between tooth and base



Isolate

- Water the model in warm water, then dry the surface
- Isolate with Isolat film
- Isolate only single-layered



Dispense

- Reposition the model
- Lock with bottom plate
- Fill dosage dispenser, because equal filling levels guarantee the optimal mixing ratio (17g powder / 11 ml liquid for a prosthesis)



Fill

- Mix Castdon Monomer and Polymer homogeneously in the mixing bowl
- Fill in a thin stream through one channel only



Polymerise

- Cure in Polymax for 25 minutes at 45°C
- Fill the water level up to the edge of the flask
- The resin-filled funnels should not be submerged






Remove

- After polymerisation, loosen the bottom plate
- Demould the model
- The denture is easily detached from the model



Fine work

- The rework is reduced to removal of surpluses, fine work with sandpaper
- ... and polishing – done!

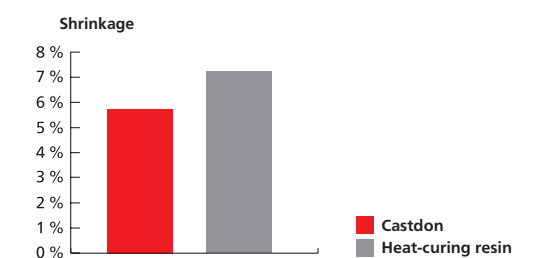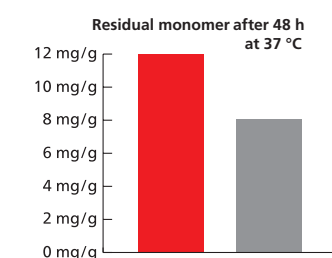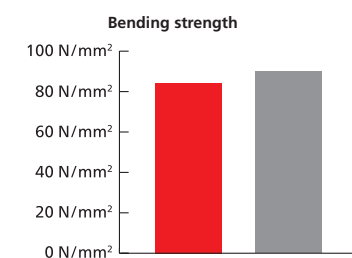

Castdon. The pouring system for full and partial dentures.

Autopolymerizing, cadmium-free PMMA-pouring resin with barbiturate curing system. Proven biocompatibility, minimum content of rest monomer, high break and color resistance. Perfectly appropriate for the production of full and partial dentures due to outstanding pouring qualities and minimum shrinkage.

- + Break-proof and color stable
- + Excellent flow characteristics
- + Lowest shrinkage

Castdon Starter Kit

Comprehensive complete package with all necessary system components for the preparation of full and partial dentures. Contains material for the preparation of about 8 full dentures. All listed items are also sold separately.

Content:

- 1 x Castdon flask (D3803)
- 1 x Castdon Monomer (500 ml) (D3865)
- 1 x Castdon Polymer (750 g) (D3873 / D3884)
- 1 x Castdon Bonding (20 ml) (D3877)
- 1 x Siliform (2 x 850 ml) (D38605)
- 1 x Isolat film (1 l) (D4304)
- 1 x Canal piercing instrument (D3872)
- 1 x Castdon dosing aids (2 pieces) (D38092)

pink transparent	D3790
pink opaque	D3792

Castdon Bonding

Liquid for the perfect bonding between highly vulcanised denture teeth and auto polymer resins.

20 ml bottle	D3877
--------------	-------

Castdon Polymer

Powder for Castdon resins based on polymethyl methacrylate copolymer.

750 g tub	transparent	D3879
	pink transparent	D3873
	pink opaque	D3884
1.2 kg tub	pink transparent	D3874
	pink opaque	D3887
4.0 kg bucket	pink transparent	D3876



The Castdon prosthesis system is perfect for repairing prostheses made with FotoDent® denture.



Castdon Monomer

Mixing liquid for Castdon resins on the basis of methyl methacrylate.

150 ml bottle	D3905
500 ml bottle	D3865
1 l bottle	D3867



Dosing aids

2 measuring cups for dosing the Castdon resin in the optimal mixing ratio.

Dosing aids	D38092
-------------	--------



Castdon flask

Professional pouring flask with transparent upper part, stainless steel bottom plate with integrated chambers for teeth dewaxing.

Caston flask incl. canal piercing instrument	D3803
----------------------------------------------	-------



Canal piercing instrument

For piercing cast canals in the dorsal area.

Canal piercing instrument, Ø 8 mm	D3872
-----------------------------------	-------



Siliform

Addition-curing, shrink-free silicone for the resin pouring technique for exact detail production and faithful reproduction of the wax set-up. Slightly oily surface film enables excellent flow characteristics. Final hardness of 14 Shore A for easy flask removal.

2 x 850 ml bottles	D38605
--------------------	--------



Gelon

High quality duplicating gel on the basis of Agar-Agar. Extra hard and stable in form. Especially developed for model casting and denture pouring technique. Fast and precise work thanks to low shrinkage and short curing times.

6.0 kg bucket	green	D38602
6.0 kg bucket	transparent	D38601



Isolat film

Heat-resistant film forming insulating agent with excellent separation effect on the basis of alginate. Excellent separation effect for smooth and pore-free surfaces. Suitable for insulation of plaster/resin as well as plaster/plaster.

1 l bottle	D4304
5 l canister	D4303N



Polymax series

Pressure polymerisation units especially developed for hot and cold curing dental resins. Adjustable temperature, 90 minutes timer, variably adjustable air pressure from 3 to 6 bar. Reliable and maintenance-free technology with automatic safety gasket system, coated aluminium pressure pot, cover from stainless steel and a sturdy housing made of polyurethane. An integrated water outlet valve facilitates the use as pressure compactor with heating switched off. Material-specific parameters of temperature and time can be monitored on the digital display and the integrated manometer.

Content:

- Pressure hose
- Water drain hose

- ⊕ Solid and maintenance-free technique
- ⊕ Automatic safety-sealing-systems
- ⊕ Housing cover made of stainless steel
- ⊕ External housing made of hard-wearing PU



Inner dimensions pressure pot:
Ø 220 mm
Height 174 mm
Volume 6.6 l

Polymax 5		up to 95 °C
red	230 V / 50–60 Hz	D3423
red	115 V / 50–60 Hz	D3423A
silver	230 V / 50–60 Hz	D34231
silver	115 V / 50–60 Hz	D34231A

Polymax 1		up to 95 °C
red	230 V / 50–60 Hz	D3429
red	115 V / 50–60 Hz	D3429A
silver	230 V / 50–60 Hz	D34291
silver	115 V / 50–60 Hz	D34291A

Polymax 1		up to 120 °C
red	230 V / 50–60 Hz	D3428
red	115 V / 50–60 Hz	D3428A
silver	230 V / 50–60 Hz	D34281
silver	115 V / 50–60 Hz	D34281A



The Polymax 5 provides enough space for 1 large articulator, 3 brass flasks with clamps or 4 castdon flasks

Polymax 3		up to 95 °C
red	230 V / 50–60 Hz	D3422
red	115 V / 50–60 Hz	D3422A
silver	230 V / 50–60 Hz	D34221
silver	115 V / 50–60 Hz	D34221A



Inner dimensions pressure pot:
Ø 160 mm
Height 110 mm
Volume 2.2 l



Inner dimensions pressure pot:
Ø 160 mm
Height 300 mm
Volume 6.0 l

Resin pouring technique

Resin pouring technique Wax boil-out units

Labormat series

Sturdy boil-out units from stainless steel for up to 6 resp. 12 flasks. Steplessly adjustable temperature from 65–95 °C and manual regulation of spray times of up to 6 minutes. A user-friendly moment switch and an integrated handspraygun add to convenient operation. For easy every-day work the unit has a 24-hour timer. The Labormat can equally be used for the curing of hot polymerizing materials.

Labormat TH

- Content:
- Cleaning brush
 - Flask supporting plate

230 V / 50 Hz	D3603N
---------------	--------



Labormat SD

- Content:
- Cleaning brush
 - Flask supporting plate
 - Water run-off / supply hose
 - Castors
 - Collection tank for waste water

230 V / 50 Hz	D3601
---------------	-------



Don't forget!

Always make sure that there is enough water in the Labormat. A too low water level leads to the activation of the safety switch.

Labormat TH	3.2 kW at 230 V / 50 Hz	65–95 °C	15 l	38 kg	440x580x360 mm
Labormat SD	3.1 kW at 230 V / 50 Hz	65–95 °C	50 l	75 kg	850x630x600 mm

Resin pouring technique Accessories

Flask basket

Made of stainless steel, with plastic handle for 3 flasks. You need 2 baskets for the Labormat TH and 4 baskets for the Labormat SD.

Flask basket	D3604
--------------	-------



Filter mats (18x60 cm)

For collecting wax residues, plaster residues etc. when boiling-out in Labormat units. You need one filter mat for the Labormat TH, two mats for the Labormat SD.

10 pieces	D36018
-----------	--------



The use of filter mats saves cleaning time

Unisol E

Descaling agent for the cleaning of frequently used boiling-out units.

2.5 kg bucket	D4383
---------------	-------



Unisol W

Concentrated wax-solving additive for a quicker and more intense cleaning of models in wax boiling-out units (e. g. Labormat).

2 x 1 l bottle	D4392
5 l canister	D4393



Base

On castors for Labormat TH with dimensions (HxWxD) 435x590x370 mm.

Base	D3607
------	-------



Let there be light.

Equipment and materials for light polymerisation.

Light-curing materials are among the most common helpers in the dental daily laboratory use. Without these and the appropriate equipment many elementary work steps cannot be adequately performed. This is why we constantly research innovations in the field of light curing. The replacement of light tubes by LEDs in our units is a shining example for simplified, more economical work processes.

Plaquit

Light-curing, one-component benzoyl peroxide-free lacquer with several application possibilities. For the coating of resin build-ups on splints, resin parts on full and partial dentures as well as temporaries in the crown and bridge technique. Polymerisation in all standard laboratory units with a spectral range of 300-400 nm, e.g. EyeEvolution®.

20 ml bottle	D38124
50 ml bottle	D38123



NanoVarnish

Transparent, light-curing one-component lacquer. Thanks to specially structured nano-particles the material creates an ideal surface energy minimising the adhesion of germs. Extremely thin layers and optimum flowing characteristics facilitate the application. Curing in light polymerisation units with the spectral area 350–500 nm, like e.g. EyeEvolution®.


20 ml bottle incl. 10 MicroBrushes	D38900
------------------------------------	--------

Tip from the pro

NanoVarnish
Especially suitable for functional areas like implants, bars, attachments.



NanoVarnish vs. conventional lacquers

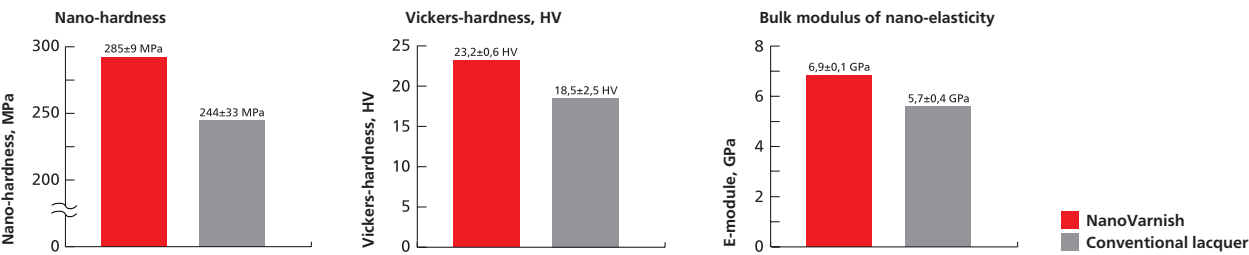

Light polymerisation Materials

Cenit L

Light-curing PMMA resin with facilitated application: Dosing aid enables applications of up to 2 mm layer thickness. For the production of space maintainers of medical and bleaching splints. Block out of undercuts or mounting of broken teeth or to fill up air bubbles in the working models. Polymerises in each standard laboratory light unit with a spectral range of 300–500 nm, e. g. EyeEvolution®.

4 application syringes (2.5 g each) incl. 8 dosing tips	D3815
---------------------------------------------------------	-------




With Cenit L, a layer thickness of up to 2 mm can be cured



Lightplast Base Plates

Pre-formed, light curing upper jaw (UJ) resin base plates for the production of custom-made impression trays, bite patterns and resin bases for the tooth set-up for complete dentures. Polymerisation in all standard laboratory units with a spectral range of 300–500 nm, e.g. EyeEvolution®.

50 pieces, pink	D38135R
50 pieces, nature	D38135N

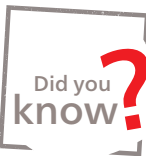

Lightplast Base Plates
Our 3D alternative: FotoDent® tray.



Individual trays, the classical application for Lightplast Base Plates



Light polymerisation Units

EyeVolution® series

The EyeVolution® light-curing units base upon the latest developments in the field of LED technique. They distinguish through their simple and intuitive handling, and due to their compact dimensions they are most suitable also for individual work stations. Most modern LED components generate exact wavelengths with extremely high light performances, thus providing for an optimum polymerisation with shortened curing times and minimum thermal strain on the material. Due to an almost unlimited durability of the LEDs the exchange of the light sources during maintenance is made redundant.

Did you know?

Different photo-initiators used in medical technology need different light sources. The light sources of the EyeVolution® series are placed according to the fields of application of these units.

- ⊕ High hardness efficiency
- ⊕ Effective deep polymerisation
- ⊕ Quiet operation
- ⊕ Low thermal stress

EyeVolution® MAX

Light polymerisation unit especially developed for all commonly available composites like in- and onlays or vestibular and occlusal veneers. Rotary reflective plate with 180° automatic swivelling and optimised LED-positioning with 6 bluelight and 1 UV light LEDs. Very short polymerisation time, e.g. thicknesses of up to 2 mm in three minutes.

EyeVolution® MAX	D38101
------------------	--------

⌚
1 / 2 / 3 / 5 / 10 min.



The reflective rotary table avoids formation of shades during polymerisation



Light polymerisation Units

EyeVolution®

Light polymerisation unit especially developed for all commercially available base materials like build-up resins of thermoformed splints, individual tray material, sealing lacquers and blockout resins. The built-in LEDs reach exactly the wavelength ranges that are needed for the activation of resins.

EyeVolution®	D38100
--------------	--------

⌚
1 / 2 / 3 / 5 / 10 min.



Positioning of the light source according to application



MultiSpot®

Polymerisation handset for lightcuring laboratory materials such as composites, splint materials, adhesives and bonding materials. Simultaneous use of UV and blue light LED with integrated performance control for all common wavelengths. A real relief in everyday laboratory work due to handy and wireless appliance.

Content:

- 1 battery charger
- 1 visor ring

⌚
0:33 min.

MultiSpot®	D38110
------------	--------



"Third hand" in dental labs: MultiSpot®



EyeVolution® MAX	100–240 V / 50–60 Hz, 0.7 A	max. 30 W	385–390 (UV light) and 465–470 nm (bluelight)	1.8 kg	205x205x255 mm
EyeVolution®	100–240 V / 50–60 Hz, 0.7 A	max. 18 W	385–390 (UV light) and 465–470 nm (bluelight)	1.5 kg	205x205x255 mm
MultiSpot®	100–240 V	2.73 W	385 (UV light) and 465 nm (bluelight)	140 g	140x30x30 mm

Dublisil®

Addition-curing vinyl polysiloxanes for all application fields in dental duplication technique, with final hardness 15/20/22/30 Shore A. Dublisil® 15 and 20 are also available as speed versions which allow deflasking after only 10 minutes. Due to the used platinum hardener the forms are virtually shrinkage-free and show an outstanding storage stability. Excellent mechanical characteristics provide for a very smooth deflasking. Negative forms can be filled several times without losing any precision. The surface of the duplication model is smooth without any additional treatment.

- ⊕ High form constancy
- ⊕ Long-term stability
- ⊕ Easy to deflask, without loss of accuracy
- ⊕ Smooth surface without additional treatment

Dublisil® 15 speed

Highspeed duplicating silicone for all kinds of duplications with a final hardness of 15 Shore A. Very easy to deflask, with highly flexible final consistency.

2 x 850 ml bottle	D4815
2 x 5.1 l canister	D4816
2 x 21.2 l canister	D4817



Dublisil® 20 speed

Universal, ultrafast duplicating silicone for all kinds of duplications with a final hardness of 20 Shore A. Especially suitable for stable duplicating forms e. g. attachment and milling techniques through the medium final hardness.

2 x 850 ml bottle	D4515
2 x 5.1 l canister	D4516
2 x 21.2 l canister	D4517



Dublisil® 15

Thinflowing standard material with a final hardness of 15 Shore A for the complete duplication technique. Very easy to deflask with highly flexible final consistency.

2 x 850 ml bottle	D4801
2 x 5.1 l canister	D4806
2 x 21.2 l canister	D4807



Dublisil® 20

Universal duplicating silicone for all kinds of duplications with a final hardness of 20 Shore A. Especially suitable for stable duplicating forms, e. g. attachment and milling techniques through the medium final hardness.

2 x 850 ml bottle	D4527
2 x 5.1 l canister	D4528
2 x 21.2 l canister	D4529



Dublisil® 22 plus

Duplication silicone with a final hardness of 22 Shore A, with a high detail accuracy and high-grade mechanical characteristics which guarantee a problem-free deflasking due to excellent tensile strength and break and tear resistance.

2 x 850 ml bottle	D4533
2 x 5.1 l canister	D4531
2 x 21.2 l canister	D4532



Dublisil® 30

Duplication silicone with a high final hardness of 30 Shore A which is characterised by its high tear resistance. Ideally applicable for highly precise embedding material e. g. telescope or attachment works. Furthermore it can be used also in flaskless duplicating technique.

2 x 850 ml bottle	D4500
2 x 5.1 l canister	D4502
2 x 21.2 l canister	D4503



Tip

from the

pro

Dublisil® speed versions
Are particularly suitable for the processing in the Dosper evo. If mixed by hand, due to the short processing time, there will only be enough time to fill one flask.

Dublisil® 15 speed	3:00 min.	10:00 min.*	> 99.8 %	15 Shore A	turquoise	medium – type 2
Dublisil® 20 speed	3:00 min.	10:00 min.*	> 99.8 %	20 Shore A	blackberry	medium – type 2
Dublisil® 15	6:00 min.	25:00 min.*	> 99.8 %	15 Shore A	green	medium – type 2
Dublisil® 20	6:00 min.	25:00 min.*	> 99.8 %	20 Shore A	pink	medium – type 2
Dublisil® 22 plus	5:00 min.	25:00 min.*	> 99.6 %	22 Shore A	red	medium – type 2
Dublisil® 30	6:00 min.	25:00 min.*	> 99.8 %	30 Shore A	purple	medium – type 2

*incl. processing time

Laboratory silicones Modelling silicones

Laborsil pro

Kneadable lab putty with a final hardness of 70 or 90 Shore A, made of addition-curing silicone. The smooth and non-sticky initial consistency enables an easy and clean mixing in a 1:1 ratio. Laborsil pro is shrinkage-free and creates a highly precise impression and exact reproduction of details. Furthermore the material is heat resistant to up to 200 °C.

Laborsil pro 70


Kneadable laboratory putty with a medium final hardness made of addition-curing silicone. Especially suitable for matrix or control key. Smooth, non sticky initial consistency, mixing in ratio of 1:1. Shrinkfree results with exact reproduction of details.

2 x 450 ml	(approx. 2 x 750 g)	D38390
2 x 1.2 l	(approx. 2 x 1.9 kg)	D38391
2 x 3.45 l	(approx. 2 x 5 kg)	D38392

Laborsil pro 90

Extremely hard, kneadable laboratory putty with smooth single components for easy kneading. Especially suitable for model material for repairs, enlargements and bite fixations. Smooth, non-sticky initial consistency, mixable in a ratio of 1:1. Shrink-free results with exact reproduction of details.

2 x 450 ml	(approx. 2 x 750 g)	D38393
2 x 1.2 l	(approx. 2 x 1.9 kg)	D38394
2 x 3.45 l	(approx. 2 x 5 kg)	D38395

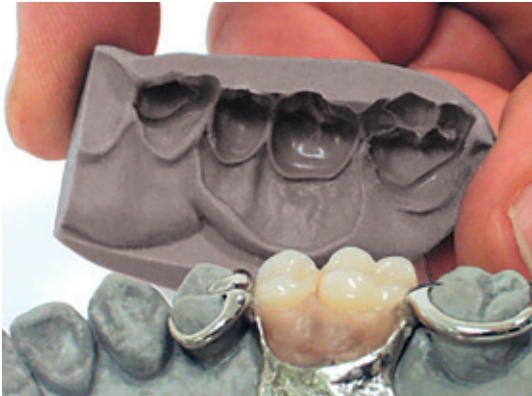

Very precise results with exact reproduction of details: Laborsil pro



Laborsil pro 70	0:30 min.	1:15 min.	6:00 min.	70 Shore A	grey	kneadable – type 0
Laborsil pro 90	0:30 min.	1:15 min.	6:00 min.	90 Shore A	grey	kneadable – type 0

Laboratory silicones Other laboratory silicones

Laborsil speed with color indicator

Versatile matrix material with color indicator. The color gradient from blue to white provides a visual control of the setting process.

2 x 50 ml double cartridges with 12 mixing tips	D38398
8 x 50 ml double cartridges	D38399

Laborsil clear

Versatile, highly transparent silicone with a final hardness of 70 Shore A. Particularly suitable as a matrix material for light-curing anterior and posterior tooth composites.

2 x 50 ml double cartridges with 12 mixing tips	D38396
8 x 50 ml double cartridges	D38397



Silfix

Precise embedding material on the basis of addition-curing silicone with perfect flow characteristics for the production of full dentures. Especially developed for the denture pressing procedure. Guarantees a safe fitting of teeth in the counter part thanks to the final hardness of over 70 Shore A.

2 x 200 ml dosing bottles with mixing beaker	D4311
4 x 50 ml double cartridges	D504310



Production of a matrix with Laborsil speed



Laborsil speed with color indicator	–	0:30 min.	2:00 min.	75 Shore A	blue to white	medium – type 2
Laborsil clear	–	0:30 min.	1:30 min.	70 Shore A	transparent	medium – type 2
Silfix	0:30 min.	2:00 min.	7:00 min.	74 Shore A	pink	low – type 3

GumQuick

Addition-curing silicone for the production of flexible gingiva masks in three versions with color adapted to natural gingiva. Exact and detailed working thanks to optimum flowing characteristics. Easy finishing with scalpel or rotating instruments. Shrinkage-free and storage stable.

GumQuick

Gingiva mask material especially developed for the C&B technique with very good tear-resistance and elongation values. Optimum flow properties for easy and safe handling.

4 x 50 ml double cartridges	D504888
-----------------------------	---------

GumQuick Set

The insulating material DuoSep included in the set also allows direct mask production in the impression.

- Content:
- 2 x 50 ml double cartridge GumQuick with 12 Mixing tips
 - 1 x 10 ml bottle DuoSep

GumQuick set	D504842
--------------	---------

Production of a gingiva mask with matrix on the model



GumQuick scan Set

Addition-curing silicone for the manufacturing of flexible and scannable gingiva masks in the CAD/CAM technique. Due to the special formulation optimal scannability or readability is achieved in commercially available dental scanners. The enclosed insulating agent DuoSep also allows a direct mask production in the impression.

- Content:
- 2 x 50 ml double cartridge GumQuick scan with 12 Mixing tips
 - 1 x 10 ml bottle DuoSep

GumQuick scan set	D504890
-------------------	---------



Developed for the digital process: GumQuick scan

GumQuick implant

Gingiva mask material especially developed for implant works. High final hardness and low flexibility for safe and exact secondary constructions. The polymerised gingiva mask shows an abrasion resistant surface.

4 x 50 ml double cartridges	D504898
-----------------------------	---------

GumQuick implant Set

The enclosed insulating agent DuoSep also allows a direct production in the impression.

- Content:
- 2 x 50 ml double cartridge GumQuick implant with 12 Mixing tips
 - 1 x 10 ml bottle DuoSep

GumQuick implant set	D504892
----------------------	---------

Production of gingiva mask in the impression with GumQuick implant



GumQuick	1:30 min.	6:30 min.	4.5 MPa	60 Shore A	low – type 3
GumQuick scan	1:30 min.	6:30 min.	4.5 MPa	65 Shore A	low – type 3
GumQuick implant	1:30 min.	6:30 min.	3.9 MPa	75 Shore A	low – type 3



Indirect denture relining with MegaBase®

MegaBase®

Permanently soft-remaining relining material on the basis of addition-curing polysiloxanes. The auto-polymerising silicone is suitable for the direct (chair-side) and for the indirect production method. The scope is very diverse and includes e. g. the relief of pressure points, protective cover after implantation or the improvement of the suction effect.

2 x 50 ml double cartridges	D5038223NR
-----------------------------	------------

MegaBase® Starter Kit

Content:

- 1 x 50 ml double cartridge MegaBase® silicone with 6 mixing tips
- 1 x 50 ml double cartridge MegaBase® 2K lacquer with 6 mixing tips
- 1 x 10 ml bottle MegaBase® Adhesive
- 3 x grinding tools
- 1 x notepad with care instructions

MegaBase® Starter Kit	D5038223
-----------------------	----------

MegaBase® 2K lacquer

Colorless, transparent, self-curing 2-component silicone lacquer for the coating of soft relining of addition-curing silicones. Very easy to use. The result: transparent-glossy and sealed surfaces.

2 x 50 ml double cartridges, 12 mixing tips	D38233
---------------------------------------------	--------

MegaBase® Adhesive

For adhesion between the addition-curing silicone and acrylate.

10 ml bottle with brush insert	D5038220NA
--------------------------------	------------



Glossy sealing of the silicone denture relining



MegaBase®	2:20 min.	20:00 min.	99.9 %	4.9 MPa	33 Shore A	medium – type 2	pink transparent
MegaBase® 2K Lack	2:15 min.	20:00 min.	–	–	25 Shore A	low – type 3	clear transparent

Odontosil

Clear-transparent, addition-curing silicone with a final hardness of 50 or 60 Shore A for the production of orthodontic elasto-positioners. The curing at room temperature permits, besides the indirect flask method, also the direct application with the double cartridge directly onto the set up model. Furthermore slight tooth corrections can be effected without any problems.

Odontosil 50

2 x 50 ml double cartridges, 12 mixing tips	D5038150
8 x 50 ml double cartridges	D5038229

Odontosil 60

2 x 50 ml double cartridges, 12 mixing tips	D5038160
8 x 50 ml double cartridges	D5038227


Odontosil 2K lacquer

Clear-transparent 2K silicone lacquer developed for coating Odontosil positioners made of addition-curing silicones. The result: transparent-glossy and sealed surfaces.

2 x 50 ml double cartridges, 12 mixing tips	D38234
---------------------------------------------	--------



Production of a silicone positioner with the Odontosil materials

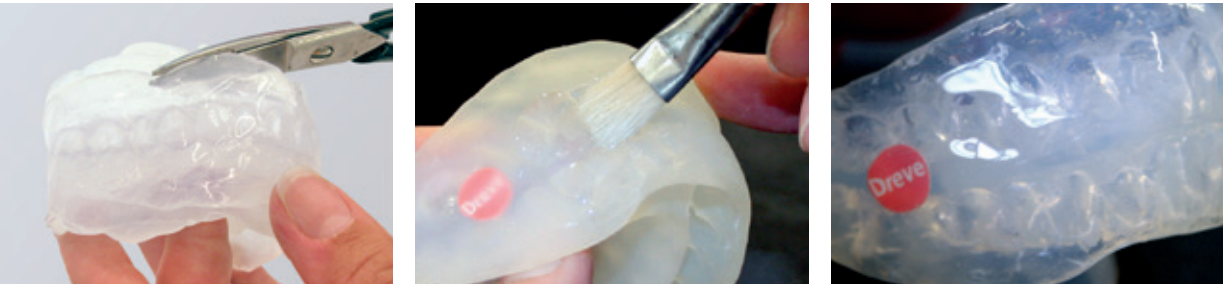

Odontosil 50	5:00 min.	30:00 min.	> 99.5 %	> 5 MPa	> 15 N/mm	< 4.0 %	50 Shore A	medium – type 2
Odontosil 60	5:00 min.	30:00 min.	> 99.5 %	> 5 MPa	> 17 N/mm	< 4.0 %	60 Shore A	medium – type 2

Laboratory silicones Mixing-/Dosing unit

Dosper evo

Automatic mixing and dosing unit for thin-flowing, addition-curing duplication silicones. Due to its innovative adjustable flow it is suitable for all duplication silicones with a mixing ratio 1:1. The result: no unwanted emptying of the reservoirs and mixing of components. The unit is constantly ready for use and very economical and is profitable also for smaller labs.

Content:

- 5 mixing tips, pink
- 2 bottle adaptors

red	230 V / 50–60 Hz	D5050
silver	230 V / 50–60 Hz	D50501
red	115 V / 50–60 Hz	D5050A
silver	115 V / 50–60 Hz	D50501A

- + Reduced formation of shades thanks to lighting of the working area
- + Shortened and safe dispensing process
- + Improvement of the material quality values
- + Low thermal strain



Dosper evo
When connecting the original containers please pierce a small hole into them in order to avoid different filling volumes.
... or simply use the material reservoirs!

Dosper evo	230 V / 115 V	250 W	12 kg	320x225x 370 mm

Laboratory silicones Accessories

Material reservoirs for Dosper M3-M5 and Dosper evo

Transparent container made of sturdy plastic for duplication silicone, equipped with quick coupling for easy installation. Filling volume approx. 1.7 l per reservoir.

2 pieces, Dosper evo	D5022
2 pieces, Dosper M3–M5	D5021



Connection hoses for Dosper M3-M5 and Dosper evo

For the safe connection between small and large DUBLISIL® canisters and the Dosper units.

2 pieces, Connection hoses for canister (5.1 l)	Dosper evo	265 mm	D49366
2 pieces, Connection hoses for canister (21.2 l)	Dosper evo	475 mm	D49367
2 pieces, Connection hoses for canister (5.1 l)	Dosper M3–M5	265 mm	D49362
2 pieces, Connection hoses for canister (21.2 l)	Dosper M3–M5	475 mm	D49363



Storage rack for Dosper M3–M5 and Dosper evo

For storage of small (5.1 l) or large (21.2 l) DUBLISIL® containers above the Dosper.

small, 371 x 295 x 214 mm (H x W x D)	D49364
big, 665 x 533 x 322 mm (H x W x D)	D49365



Clean duplication process with Dosper evo



- ⊕ Safe fitting of the negative due to retention part
- ⊕ Fixture with conical cuff
- ⊕ Top, cuff and base part are made of durable Makrolon plastic



Height adjustable flask

Dublisil® flask

Height adjustable flask for the economic silicone duplication technique for models with or without base. Highly economical through variably adjustable interior part and therefore best adjustment to the corresponding model situation.

Dublisil flask, incl. foam spacer	D3805
-----------------------------------	-------




Economic duplications thanks to Dublisil® flask

Foam spacer

For Dublisil® flasks, encloses the model base for more economical duplications.

25 pieces	D38076
-----------	--------

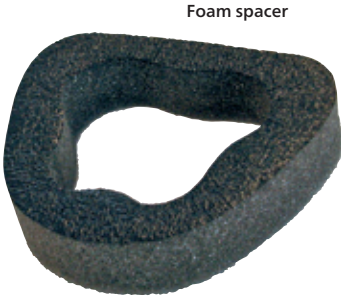

Foam spacer

DuoSep

Insulation agent between addition-curing silicone and polyether materials.

10 ml bottle with brush insert	D4306
--------------------------------	-------



Mixing tips

Mixing ratio	1:1
Inner-Ø	5.4 mm
Color	pink



Dosper evo / M5
Silfix

40 pieces	D4970X1
100 pieces	D49701X1

Mixing ratio	1:1
Inner-Ø	6.5 mm
Color	green

Regofix
Laborsil clear
Laborsil speed
MegaBase®
Odontosil

40 pieces	D4970X2
100 pieces	D49701X2



Mixing ratio	1:1
Inner-Ø	5.4 mm
Color	transparent



Dosper M4

40 pieces	D4960
100 pieces	D49601

Mixing ratio	1:1
Inner-Ø	3.2 mm
Color	light blue

GumQuick
GumQuick scan
GumQuick implant
MegaBase® 2K Lack
Odontosil 2K Lack

40 pieces	D4970X3
100 pieces	D49701X3



DS 50 Injector 1:1

For applications of dental laboratory silicones from S 50 double cartridges.

Mixing ratio	1:1
DS 50 Injector 1:1	D5007



More products and accessories

Modralit® 3K Set

Castable model stump material based on polyurethane for the production of precision models. Due to its edge strength, Modralit® 3K is excellently suited for the precious metal and base metal technology and is particularly suitable for step and partial crown preparations. The set is enough for about 10 complete tooth rims, which can then be coated with commercial painting systems.


Content:

- 2 x 100 ml bottles base
- 1 x 100 ml bottle hardener
- 1 x 400 g tub filler material

Modralit® 3K Set	D3926
------------------	-------

Refill packs for Modralit® 3K Sets

100 ml	Base	D39261
100 ml	Hardener	D39262
400 g	Filler	D39263

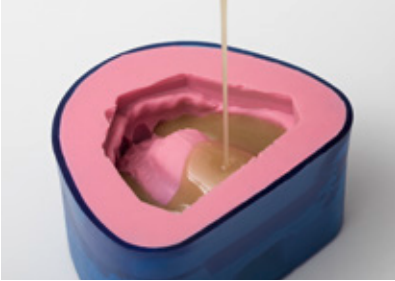

The precise alternative to the plaster model:
Modralit® 3K



Separator

Agent for the separation of non-compatible materials like Modralit® 3K from polyether impressions, A-silicones from C-silicones and identic materials, for example A-silicones from A-silicones.

50 ml bottle with brush insert	D4305
--------------------------------	-------

GipsEx

Chemical plaster solvent for use e. g. in the resin technique. GipsEx is especially effective and consequently time-saving when cleaning resin dentures, mixing beakers etc. Suitable for use in ultra-sonic units.

2 x 1 l bottle	D4395
1 x 5 l canister	D4396



Modralit® 3K	~ 2:00 min	~ 0.07 %	ivory

More products and accessories

UltraTrimm

Flexible grinding wheels for the rough to fine elaboration of EVA thermoforming blanks for sports mouthguards (e. g. Drufosoft®), soft-remaining relining silicone (e. g. MegaBase®), positioner silicone (e. g. Odontosil) as well as hard and soft/hard thermoforming blanks at a maximum of 15.000 revolutions / minute.

UltraTrimm set


Set with several grinding wheels in different selectable hardnesses. Each set additionally contains 4 support discs as well as 1 mandrel.

UltraTrimm set light, 12 grinding wheels	D43753
UltraTrimm set medium, 12 grinding wheels	D43752
UltraTrimm set heavy, 6 grinding wheels	D43751

UltraTrimm economy pack

Storage pack with 80 resp. 40 grinding wheels in different selectable hardnesses.

UltraTrimm Set light, 80 grinding wheels	D43763
UltraTrimm Set medium, 80 grinding wheels	D43762
UltraTrimm Set heavy, 40 grinding wheels	D43761

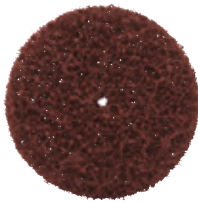

UltraTrimm light
(Ø 28 mm / thickness 7 mm)



UltraTrimm medium
(Ø 28 mm / thickness 7 mm)



UltraTrimm heavy
(Ø 28 mm / thickness 14 mm)



Mandrel and
support disc

Grinding cap holders and grinding caps

Rotating, economic and interchangeable corundum grinding caps with a grain of 150 µm for the elaboration of EVA thermoforming blanks for sports mouthguards (e. g. Drufosoft®), soft-remaining relining silicone (e. g. MegaBase®) as well as positioner silicone (e. g. Odontosil) at a maximum of 20,000 revolutions / minute.

Size	1 x Grinding cap holder			5 x grinding caps	
1	round; Ø 7 mm D38230T1		◀ compatible with ▶	round; Ø 8.5 mm D38230S1	
	pointed; Ø 5 mm D4458			pointed; Ø 6.5 mm D4461	
2	round; Ø 10 mm D38230T2		◀ compatible with ▶	round; Ø11.5 mm D38230S2	
	pointed; Ø 7 mm D4459			pointed; Ø 8.5 mm D4462	
3	round; Ø 13 mm D38230T3		◀ compatible with ▶	round; Ø 14.5 mm D38230S3	

More products and accessories

Mixing tips

Mixing ratio 1:1
Inner-Ø 5.4 mm
Color pink

Dosper evo / M5
StoneBite®
StoneBite® scan
Silfix

40 pieces	D4970X1
100 pieces	D49701X1



Mixing ratio 1:1
Inner-Ø 4.2 mm
Color yellow

Dynax® light


40 pieces	D4970X5
100 pieces	D49701X5



Mixing ratio 4:1/10:1
Inner-Ø 3.2 mm
Color light blue

Fixtemp C&B 4:1
Fixtemp C&B 10:1

40 pieces	D4970X4
100 pieces	D49701X4

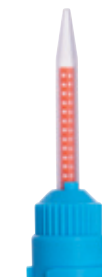

Mixing ratio 1:1
Inner-Ø 3.2 mm
Color light blue

GumQuick
GumQuick scan
GumQuick implant
MegaBase® 2K Lack
Odontosil 2K Lack

40 pieces	D4970X3
100 pieces	D49701X3



Mixing ratio 1:1
Inner-Ø 6.5 mm
Color green

Dynax® clear
Dynax® heavy body
Dynax® mono
Regofix
Laborsil clear
Laborsil speed
Megabase®
Odontosil

40 pieces	D4970X2
100 pieces	D49701X2



Mixing ratio 1:1
Inner-Ø 2.5 mm
Color brown

for double-chamber syringes

40 pieces	D4970X6
100 pieces	D49701X6



Mixing ratio 1:1
Inner-Ø 5.4 mm
Color transparent

Dosper M4

40 pieces	D4960
100 pieces	D49601



Mixing ratio 1:1
Inner-Ø 4.2 mm
Color mint green

Dynax® light

40 pieces	D4970X7
100 pieces	D49701X7



More products and accessories

Tips

For mixing tips for intra-oral application.

transparent (e. g. for Dynax® clear)
yellow (e. g. for Dynax® light)

40 pieces, transparent	D49602
40 pieces, yellow	D49604



Application tips

For mixing tips for intra-oral application. Flattened form for bite registration materials.

transparent (e. g. for StoneBite®)

40 pieces, transparent	D49603
------------------------	--------



Bayonet-ring

Attachment for the use of 380 ml cartridges.

Bayonet-ring, yellow	D4980
----------------------	-------



Application tips

Top pieces for dosing syringes.

Inner-Ø 0.9 mm
black (e. g. for Cenit L)

15 pieces, black	D2308
------------------	-------



Mixing tip

Mixing ratio 5:1
Inner-Ø 6.0 mm
Color yellow

Zerosil® soft
Dynax® heavy body
Dynax® mono
Dynax® putty

15 pieces	D4982
50 pieces	D4981



The colors of the closure caps on the double cartridges indicate the mixing tips that belong to the respective product.

DS 50 Injector 1:1

For applications of dental laboratory silicones from S 50 double cartridges.

Mixing ratio 1:1	
DS 50 Injector 1:1	D5007

DS 50 Injector 4:1/10:1

For the application of C&B materials from S 50 double cartridges.

Mixing ratio 4:1/10:1	
DS 50 Injector 4:1/10:1	D5008



Index

#		F		N	
3D Protector inlays	29	Fillin	35	NanoVarnish	46
A		Filter mats	45	O	
Application tips – black	65	Finishing Liquid	35	Odontosil 2K lacquer	57
Application tips – transparent	14, 65	Fixtemp C&B 4:1	15	Odontosil 50	57
B		Fixtemp C&B 10:1	15	Odontosil 60	57
Bandelin RK 100 H	20	Flask basket	45	P	
Base	45	Foam spacer	60	PCU LED N ₂	21
Bayonet-ring	14, 65	Foil scissors	34	Piercing cutter	34
Beaker	20	FotoClean	21	Plaquit	46
BioBleach® hard	24	FotoDent® biobased model	19	Polymax 1 (95 °C / 120 °C)	43
BioBleach® soft	24	FotoDent® cast	18	Polymax 3	43
Biolon	24	FotoDent® denture	18	Polymax 5	42
Biolon for Vacformat V9 2.0	25	FotoDent® gingiva	19	Pressure polymerisation insert	37
Blank trays	25	FotoDent® guide	18	Protector adhesive	29
C		FotoDent® IBT	18	R	
Canal piercing instrument	41	FotoDent® model2	19	Regofix opak	13
Castdon Bonding	40	FotoDent® setup	18	Regofix transparent	13
Castdon flask	41	FotoDent® tray	18	S	
Castdon Monomer	41	Funky tool	37	Sample selection	
Castdon Polymer	40	G		– thermoforming blanks	24
Castdon Starter Kit	40	Gelon	41	Separator	62
Cenit L	47	GipsEx	62	Silicone grease	35
Connection hoses		Granules	36	Silfix	53
– for Dosper M3-M5/evo	59	Grinding caps	63	Siliform	41
D		Grinding cap holders	63	Stainless steel perforated lid	20
Dosing aids	41	GumQuick	54	StoneBite®	13
Dosper evo	58	GumQuick implant	55	StoneBite® scan	13
Dreve Mouthguard Label	30	GumQuick scan Set	54	Storage racks	59
Dreve Mouthguard Box	30	I		Square cutter	34
Dreve Mouthguard professional 3D		Insert basket made of stainless steel	20	T	
– Starter Set	29	Isolat film	41	Tips – transparent	14, 65
Drufolen H	25	K		Tips – yellow	14, 65
Drufolen W	25	Kombiplast	27	U	
Drufomat scan	32	L		UltraTrimm economy pack	63
Drufosmart	33	Labormat SD	44	UltraTrimm Set	63
Drufosmart scan	33	Labormat TH	44	Unisol E	45
Drufoplast	27	Laborsil clear	53	Unisol W	45
Drufoplast incl. tray handle	27	Laborsil pro 70	52	V	
Drufosoft®	26	Laborsil pro 90	52	Vacformat V9 2.0	31
Drufosoft® pro	27	Laborsil speed	53	W	
DS 50 Injector 1:1	14, 61, 65	Lightplast base plates	47	Workout set	34
DS 50 Injector 4:1 / 10:1	15, 65	M		X	
Dublisil® 15	50	Material reservoirs		X-ray opaque stripes	30
Dublisil® 15 speed	50	– for Dosper M3–M5/evo	59	Z	
Dublisil® 20	51	MaxiFresh™ cleaning tablets	36	Zerosil® soft	12
Dublisil® 20 speed	50	MegaBase®	56	Zerosil® soft	
Dublisil® 22 plus	51	MegaBase® 2K lacquer	56	– with color indicator	12
Dublisil® 30	51	MegaBase® Adhesive	56		
Dublisil® flask	60	MegaBase® Starter Kit	56		
DuoSep	61	Mixing tips – 4:1 / 10:1, light blue	15, 64		
Dynax® clear	11	Mixing tips – 5:1, yellow	14, 65		
Dynax® heavy body	10	Mixing tips – brown	64		
Dynax® light	10	Mixing tips – yellow	14, 64		
Dynax® mono	11	Mixing tips – green	14, 61, 64		
Dynax® putty	10	Mixing tips – light blue	61, 64		
E		Mixing tips – light blue	61, 64		
EyeVolution®	49	Mixing tips – mint green	14, 64		
EyeVolution® MAX	48	Mixing tips – pink	14, 61, 64		
		Mixing tips – transparent	61, 64		
		Modralit® 3K Set	62		
		Modralit® 3K Set			
		– Refill packs	62		
		MultiSpot®	49		

Medical devices

Our medical devices are registered in many countries.
Please contact us for further information.

We cannot be held responsible for false statements in this catalogue and their consequences. Delivery of our products is effected according to our general terms and conditions. Technical changes reserved. We do not take any responsibility for the authenticity of colors in this catalogue.



Dreve Dentamid GmbH
Max-Planck-Straße 31
59423 Unna/Germany

Phone +49 2303 8807-37
Fax +49 2303 8807-55
Email dentamid@dreve.de
Internet www.dentamid.dreve.de

Application technology

Phone +49 2303 8807-3404

Sales Germany

Phone +49 2303 8807-40
Fax +49 2303 8807-55

Sales Export

Phone +49 2303 8807-37
Fax +49 2303 8807-55

Design & product photos

Dreve ProDiMed GmbH
Advertising department
Einsteinstraße 38
59423 Unna/Germany



Dreve Dentamid GmbH

Max-Planck-Straße 31 • 59423 Unna/Germany
Telefon +49 2303 8807-37 • Telefax +49 2303 8807-55
dentamid@dreve.de • www.dentamid.dreve.de